



中国认可
国际互认
检测
TESTING
CNAS L1842



佛山中纺联检验技术服务有限公司
CNTAC Testing Service Co.,Ltd.(Foshan)

Testing Report	Security website: www.fcl-sz.org.cn Security code: 8654220434
Report No: ZFLJ2664271	Page 1 of 4



Applicant Information

Applicant Name : AccuMed Biotech LLC
Applicant Address : 16727 Park Row Houston TX 77084 USA
Manufacturer : Gaomi Hancheng Labor Protection Appliance Co.,Ltd

Sample Information

Sample Description : White KN95 Mask With Earloop
Brand : AccuMed
Sample Quantity : 20 pieces
Material Components : 6 layers 1.45gsm nonwoven fabric 2.25gsm meltblown 3.25gsm nonwoven fabric 4.Filled filter cotton 5.25gsm meltblown 6.45gsm nonwoven fabric

- Sample Receiving/Testing date : 2020-12-15
- Report Date : 2020-12-21
- The original sample is stucked on the last paper.

Test Performed

Judgement according to:
GB 2626-2019 Respiratory protection— Non-powered air-purifying particle respirator .
- Selected test(s) as requested by applicants. For details, refer to attached page(s).

Pronounce

The results shown in this report refer only to sample(s) tested unless otherwise stated.
All the items are performed in the standard conditions, except the noted cases.
Except for the requirement of the client, the test results and the conformity judgement of this report do not take the uncertainty of the test results into account.

Signed for and on behalf of
CNTAC Testing Service Co.,Ltd.(Foshan)

Approved by

张志荣





中国认可
国际互认
检测
TESTING
CNAS L1842



佛山中纺联检验技术服务有限公司
CNTAC Testing Service Co.,Ltd.(Foshan)

Testing Report	Security website: www.fcl-sz.org.cn Security code: 8654220434
Report No: ZFLJ2664271	Page 2 of 4



CONCLUSION

Test Item
Filtration Efficiency

Judgement
Pass



Testing Report	Security website: www.fcl-sz.org.cn Security code: 8654220434
Report No: ZFLJ2664271	Page 3 of 4



--- Test Result ---


Test Result	Requirements	Judgement
1. Filtration Efficiency <u>GB 2626-2019 Section 6.3</u>		
Air flow:85L/min,Aerosol: NaCl		
Unit: <%>		
Before pretreatment	(KN95) ≥ 95.0	Pass
sample 1#	99.6	
sample 2#	99.6	
sample 3#	99.5	
sample 4#	99.5	
sample 5#	99.4	
sample 6#	99.6	
sample 7#	99.5	
sample 8#	99.6	
sample 9#	99.5	
sample 10#	99.5	
sample 11#	99.4	
sample 12#	99.7	
sample 13#	99.6	
sample 14#	99.6	
sample 15#	99.5	
After pretreatment		
sample 1#	99.5	
sample 2#	99.4	
sample 3#	99.6	
sample 4#	99.7	
sample 5#	99.5	
Minimum	99.4	



Testing Report	Security website: www.fcl-sz.org.cn Security code: 8654220434
Report No: ZFLJ2664271	Page 4 of 4



Original Sample

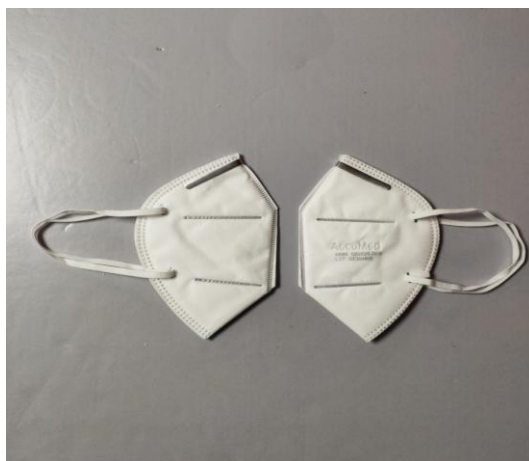

===== End of Report =====

