



中国认可  
国际互认  
检测  
TESTING  
CNAS L1842



佛山中纺联检验技术服务有限公司  
CNTAC Testing Service Co.,Ltd.(Foshan)

<b>检测报告(Testing Report)</b>	防伪查询网址(Security website): <a href="http://www.fcl-sz.org.cn">www.fcl-sz.org.cn</a> 防伪码(Security code): 1058016116
报告编号(Report No): ZFLJ2625681	Page 1 of 4



### Applicant Information (客户信息)

Applicant Name(委托单位) : AccuMed Biotech LLC  
Applicant Address(委托单位地址) : 16727 Park Row Houston TX 77084 USA  
Manufacturer(生产单位) : Gaomi Hancheng Labor Protection Appliance Co.,Ltd

### Sample Information (样品信息)

Sample Description(样品描述) : KN95 Mask  
Sample Quantity(样品数量) : 49 pieces (49 个)  
Sample Grade(样品等级) : BFE95  
Colour(颜色) : white  
Material Components(原料成分): Non-woven fabric melt-blown fabric hot air cotton (无纺布 熔喷布 热风棉)

- Sample Receiving date (收到样品的日期) : 2020-05-28  
- Report date (报告日期) : 2020-06-02  
- The original sample is stucked on the last paper. (送检样品原样粘贴在本报告的末页)

### Test Performed (检验标准)

Judgement according to(评定依据):  
GB 2626-2006 Respiratory protective equipment-Non-powered air-purifying particle respirator  
(呼吸防护用品 自吸过滤式防颗粒物呼吸器)

- Selected test(s) as requested by applicants. (本实验室根据客户要求完成以下检测内容)

### Pronounce (声明)

The results shown in this report refer only to sample(s) tested unless otherwise stated.All the items are performed in the standard conditions, except the noted cases.(除非特别说明,测试结果只对来样负责。所有测试项目均在标准规定的环境下进行,有注明除外。)

Except for the requirement of the client, the test results and the conformity judgement of this report do not take the uncertainty of the test results into account.(除非客户要求,本报告检测结果及符合性判定不考虑测量结果的不确定度。)

Signed for and on behalf of

**CNTAC Testing Service Co.,Ltd.(Foshan)**

Approved by  
(批准)

张志荣





中国认可  
国际互认  
检测  
TESTING  
CNAS L1842




<b>检测报告(Testing Report)</b>	防伪查询网址(Security website): <a href="http://www.fcl-sz.org.cn">www.fcl-sz.org.cn</a>
	防伪码(Security code): 1058016116
报告编号(Report No): ZFLJ2625681	Page 2 of 4



--- Test Result (测试结果) ---

Test Result(测试结果)	Requirements(评定条件)	Judgement(判定)
<b>1. Head harness(头带) <u>GB 2626-2006 Section 6.11(条款 6.11)</u></b>		
For disposable facepiece, tension 10N for 10s(随弃式面罩, 拉力 10N, 持续 10s)		
Before pretreatment (未经处理)	No slippage, fracture shall occur(不应出现滑脱或断裂)	Pass(合格)
No slippage or fracture (未出现滑脱和断裂)		
After pretreatment(经预处理)		
No slippage or fracture (未出现滑脱和断裂)		
<b>2. Filtration Efficiency(过滤效率) <u>GB 2626-2006 Section 6.3(条款 6.3)</u></b>		
Air flow: 85L/min, Aerosol: NaCl(气体流量 85L/min, 气溶胶颗粒: NaCl)		
Unit(单位): <%>		
Before pretreatment(未经处理)	(KN95) ≥ 95.0	Pass(合格)
sample 1#(样品 1)	99.5	
sample 2#(样品 2)	99.5	
sample 3#(样品 3)	99.3	
sample 4#(样品 4)	99.4	
sample 5#(样品 5)	99.4	
sample 6#(样品 6)	99.5	
sample 7#(样品 7)	99.5	
sample 8#(样品 8)	99.7	
sample 9#(样品 9)	99.6	
sample 10#(样品 10)	99.5	
After pretreatment(经预处理)		
sample 1#(样品 1)	99.5	
sample 2#(样品 2)	99.4	
sample 3#(样品 3)	99.5	
sample 4#(样品 4)	99.4	
sample 5#(样品 5)	99.4	
<b>3. Exhalation Resistance(呼气阻力) <u>GB 2626-2006 Section 6.6(条款 6.6)</u></b>		
Air flow: 85 L/min(通气量 85 L/min)		
Unit(单位): <Pa>		
Before pretreatment(未经处理)	≤ 250	Pass(合格)
sample 1#(样品 1)	89.0	
sample 2#(样品 2)	92.3	
After pretreatment(经预处理)		
sample 1#(样品 1)	100.7	
sample 2#(样品 2)	109.0	

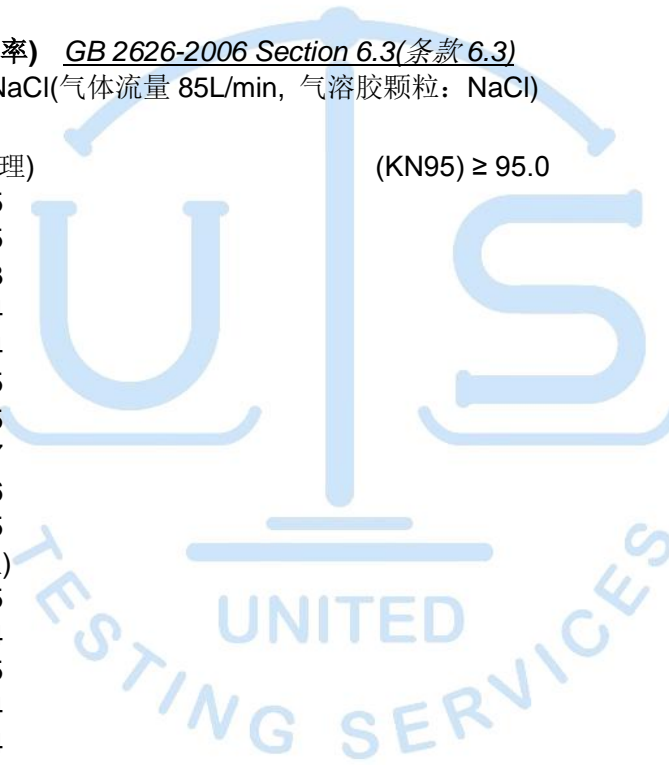




中国认可  
国际互认  
检测  
TESTING  
CNAS L1842



佛山中纺联检验技术服务有限公司  
CNTAC Testing Service Co.,Ltd.(Foshan)

<b>检测报告(Testing Report)</b>	防伪查询网址(Security website): <a href="http://www.fcl-sz.org.cn">www.fcl-sz.org.cn</a>	
	防伪码(Security code): 1058016116	
报告编号(Report No): ZFLJ2625681	Page 3 of 4	

Test Result(测试结果)

Requirements(评定条件)

Judgement(判定)

**4. Inhalation Resistance(吸气阻力)** GB 2626-2006 Section 6.5(条款 6.5)

Air flow: 85 L/min(通气量 85 L/min)

Unit(单位): <Pa>

Before pretreatment(未经处理)

≤ 350

Pass(合格)

sample 1#(样品 1) 101.5

sample 2#(样品 2) 110.7

After pretreatment(经预处理)

sample 1#(样品 1) 108.9

sample 2#(样品 2) 115.2

Remark (备注): The filtration efficiency, inhalation resistance and exhalation resistance items are unauthorized for CNAS.  
(CNAS 未授权过滤效率、吸气阻力和呼气阻力项目。)



<b>检测报告(Testing Report)</b>	防伪查询网址(Security website): <a href="http://www.fcl-sz.org.cn">www.fcl-sz.org.cn</a> 防伪码(Security code): 1058016116
报告编号(Report No): ZFLJ2625681	Page 4 of 4



Original (送检样品原样)



==== End of Report (检测报告结束) =====

